


Intramuscular and intradermal vaccination trainer

PRODUCT NUMBER: 2514

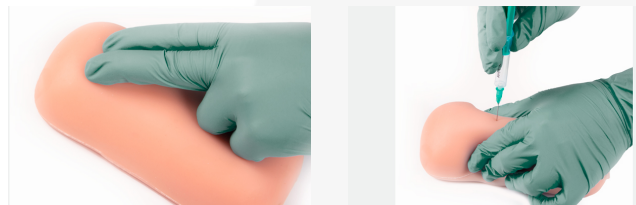

DESCRIPTION:

This trainer enables **intramuscular injection** to be practiced, as is typically done when **vaccinating**. The model can be used as a vaccination trainer, but also as a normal IM injection trainer and as a subcutaneous injection trainer.

SKILLS:

- Intramuscular injection
- Subcutaneous injection
- Syringe and needle handling

MORE THAN SIMULATORS S.L.

Calle de la Hoya, 4 2º
44570 Calanda (Teruel)
Spain

CONTACTO

sales@morethansimulators.com
(+34) 936 06 77 05 / (+351) 21 123 85 18
<https://morethansimulators.com>



CHARACTERISTICS:

GENERAL:

- Realistic anatomical landmarks
- Dimensions: 17.5 x 8.5 x 5.0 cm
- Weight: 300 g

FEATURES:

- The trainer can be buckled and then simulates a human upper arm. A simulated and easily palpable shoulder anatomy serves as a landmark for the puncture.
- The back of the trainer is puncture-proof so that there is no danger to the patient.
- The injection can be done with normal needles and syringes, the injection of liquid is possible.
- The model can be dismantled so that the liquid can easily be removed after the exercise.
- Thanks to the new material, the arm feels lifelike and the puncture points are not visible, so that many hundreds of punctures are possible before the vaccination trainer has to be replaced.

INCLUDE:

- 1 Vaccination trainer
- 1 Instructions manual

MORE THAN SIMULATORS S.L.

Calle de la Hoya, 4 2º
44570 Calanda (Teruel)
Spain

CONTACTO

sales@morethansimulators.com
(+34) 936 06 77 05 / (+351) 21 123 85 18
<https://morethansimulators.com>